

## SKIN GRAFTING PRODUCTS ORDERING INFORMATION

### Skin Graft Mesher

Prod. No.	Description
7701-00	Zimmer Skin Graft Mesher <sup>1</sup>
7701-01	Zimmer Skin Graft Mesher Complete <sup>2</sup>
7701-02	Zimmer Skin Graft Mesher w/Two Cutters <sup>3</sup>
7701-03	Zimmer Skin Graft Mesher Autoclave Case
7701-22	Zimmer Skin Graft Mesher Ratchet Handle

### Skin Graft Mesher Cutters

Prod. No.	Description
7703-15	1.5:1 Expansion Ratio
7703-20	2:1 Expansion Ratio
7703-30	3:1 Expansion Ratio
7703-40	4:1 Expansion Ratio

### Skin Graft Mesher Carriers\*

Prod. No.	Description
7708-00-10	<i>Dermacarriers</i> Length 8" (20.3cm)
7716-00-10	<i>Dermacarriers</i> Length 16" (40.6cm)

\* *Dermacarriers* Skin Graft Carriers (7708-00-10 and 7716-00-10) may be used only with the Zimmer Skin Graft Mesher.

### Companion Products

Prod. No.	Description
8801-00	Zimmer Air Dermatome Complete <sup>4</sup>
8801-01	Zimmer Air Dermatome Handpiece <sup>5</sup>
8821-00	Zimmer Electric Dermatome Complete <sup>6</sup>
8821-01	Zimmer Electric Dermatome Handpiece <sup>7</sup>
8800-00-10	Zimmer Dermatome Blades (10/box)
8821-06	Power Supply <sup>8</sup>

### MeshGraft II Tissue Expansion System

Prod. No.	Description
2195-00	<i>MeshGraft II</i> Complete <sup>9</sup>
2195-01	<i>MeshGraft II</i> Instrument only

### MeshGraft II Carriers\*

Prod. No.	Description
2195-12	1.5:1 <i>Dermacarriers II</i> Skin Graft Carrier
2195-13	3:1 <i>Dermacarriers II</i> Skin Graft Carrier
2195-14	6:1 <i>Dermacarriers II</i> Skin Graft Carrier
2195-15	9:1 <i>Dermacarriers II</i> Skin Graft Carrier

\* *Dermacarriers II* Skin Graft Carriers (2195-12 through 2195-15) may be used only with the *MeshGraft II* Tissue Expansion System.


#### Footnotes

- 1 Includes mesher, autoclave case, and ratchet handle.
- 2 Includes mesher, autoclave case, ratchet handle, and cutters (1.5:1, 2:1, 3:1, 4:1).
- 3 Includes mesher autoclave case, ratchet handle, and two cutters (1.5:1, 3:1).
- 4 Includes handpiece, hose, autoclave case, width plates, and screwdriver.
- 5 Includes handpiece, hose, width plates, and screwdriver.
- 6 Includes handpiece, power supply, power cord, autoclave case, width plates, and screwdriver.
- 7 Includes handpiece, power cord, width plates, and screwdriver.
- 8 The power supply is shipped at 120 VAC. Additional fuses, 1.5A for use at 100 VAC and 0.8A for use at 220-240 VAC, are packaged with the power supply for international use.
- 9 Includes mesher, case, ratchet handle, Allen wrench, and manual.

For more information, contact your Zimmer representative, or visit us at [www.zimmer.com](http://www.zimmer.com).



## ZIMMER® SKIN GRAFT MESHER

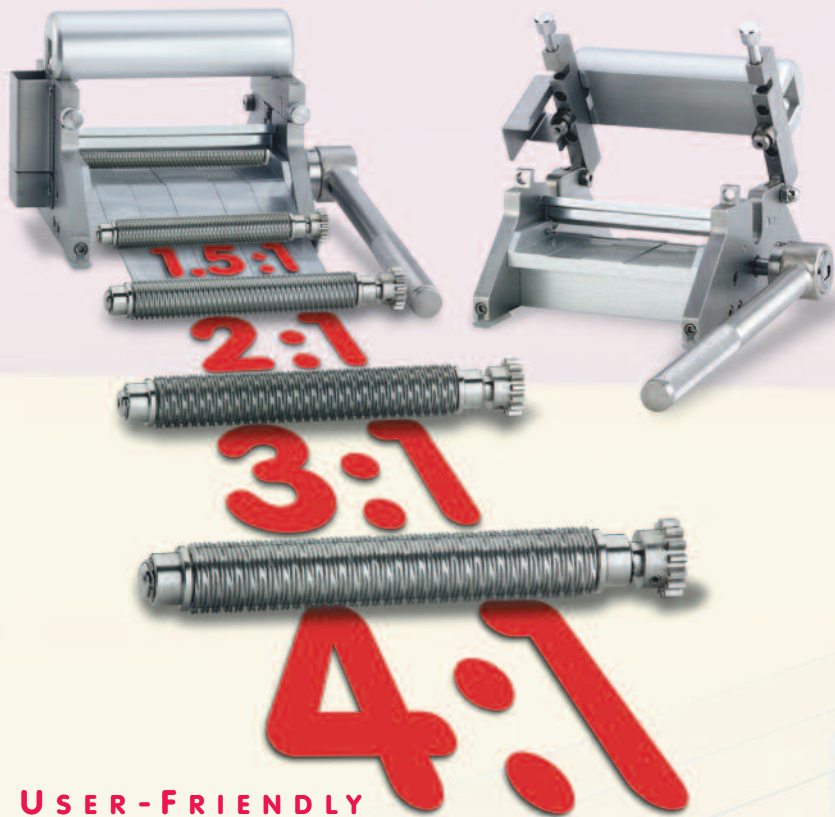

### USER-FRIENDLY

The *Zimmer* Skin Graft Mesher is designed for simple operation.

- Durable, stainless-steel construction for years of consistent performance.
- Simple to change-out cutter options.
- Easy Ratchet handle motion for smooth transport of graft carrier through the mesher.

### EFFICIENT

The Skin Graft Mesher's interchangeable cutters provide for true efficiency.

- Four common expansion ratios help to provide patient specific results by promoting proper wound drainage and epithelialization for proper healing.
- Intraoperative efficiencies are realized, without the hassle of multiple machines.

### DURABLE

The *Zimmer* Skin Graft Mesher is engineered to provide durability, convenience, and dependability in surgical skin grafting procedures:

- Convenient hinge design allows the cutters to be changed quickly or removed for easier use and cleaning.
- Built-in comb minimizes the chances of grafts sticking to and wrapping around the cutters.
- Stainless-steel case makes sterilization and transportation simple.
- Lubeless gears help to ensure a smooth performance, time after time.

### ZIMMER SKIN GRAFT MESHER LARGER DISPOSABLE CARRIERS FOR EASIER GRAFT HANDLING

Disposable carriers, available in 8" x 4" and 16" x 4", guide large, multiple, or small grafts through the mesher. The larger carriers also help to minimize the risk of damage as fragile grafts are transported to the site.

### MESHGRAFT® II TISSUE EXPANSION SYSTEM



### MESHGRAFT II SYSTEM

The *Meshgraft* II System offers four skin ratios, 1.5:1, 3:1, 6:1, and 9:1 through the use of the *Dermacarriers*® II Skin Graft Carriers.

The *Meshgraft* II instrument consists of a ratchet handle, continuous feed roller, a cutting roller, and a guidance plateau. The guidance plateau helps ensure a proper alignment of the *Dermacarriers* II Skin Graft Carriers and cutting rollers, and the ratchet handle facilitates proper advancement of the carrier.



## SKIN GRAFT MESHER

User-Friendly  
Efficient  
Durable