




**Zimmer®
Coonrad/Morrey
Total Elbow**



Interchangeability, Anterior Flange, Clinical Success

The World's Number One Selling Total Elbow¹

For complete patient matching flexibility, the Coonrad/Morrey Total Elbow is unsurpassed in interchangeability and design and is the total elbow of choice for its ease of use. Backed by more than 20 years of clinical success,^{2,4} the Coonrad/Morrey Total Elbow is indicated for primary and revision intervention, post traumatic² and rheumatoid arthritis,^{3,6} and trauma reconstruction.^{4,5} The bottom line—results. That's what ease-of-operation, enhanced fixation, and simple, precise instrumentation can provide . . . that's the Coonrad/Morrey Total Elbow.

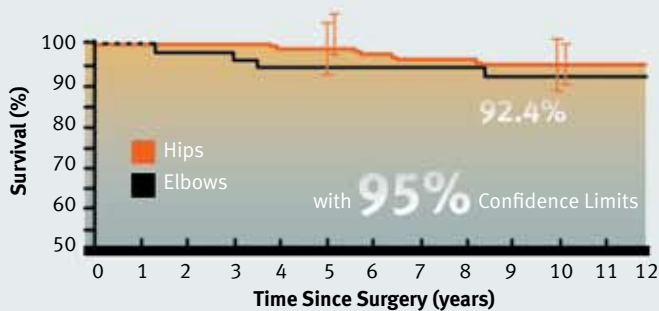

Clinical Success

In joint replacement surgery, total hip surgery offers a standard for clinical success. In comparison the Coonrad/Morrey Total Elbow has performed equally or better^{2,3} and is backed by more than 20 years of clinical success.^{2,4}





Trauma Reconstruction

The longer flange option addresses clinical situations involving distal/humeral bone loss without compromising elbow joint performance



Complete Interchangeability


Twelve Humeral Stems – 4”, 6”, and 8” lengths in extra-small, small, and regular sizes

Twelve Ulnar Stems – regular and long lengths in extra-small, small, and regular sizes

Provides comprehensive patient matching and interoperative flexibility

Two-piece locking pin simply snaps together and provides accuracy and efficiency

Easy disassembly of the locking pin enhances interchangeability in revision cases

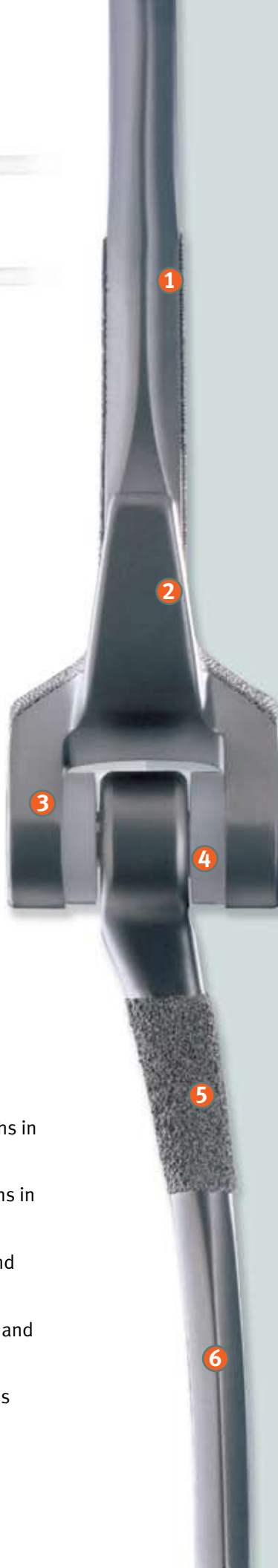

1

Triangular to match canal landscape maximizing stability by minimizing rotation

2

Proximal humeral flange

permits insertion of bone graft anteriorly to enhance thickening of bone stock where maximum stress occurs, reducing posterior displacement and torsional forces and helping protect bone/cement/implant interface

3

Titanium® Ti-6Al-4V Alloy

for durable, lightweight construction

4

Zimmer Certified Poly Bushings

with ultra-high molecular weight polyethylene, prevents metal-to-metal contact

5

Plasma spray coating

provides enhanced cement fixation; implant intended to be used with bone cement for both immediate and long-term fixation

6

Quadrangular to match canal landscape

to help minimize intramedullary rotation in the ulna

Implants

Humeral Implants

Prod. No.	Description
32-8105-025-04	Humeral Component, 4", Small
32-8105-025-06	Humeral Component, 6", Small
32-8105-025-08	Humeral Component, 8", Small
32-8105-026-04	Humeral Component, 4", Regular
32-8105-026-06	Humeral Component, 6", Regular
32-8105-026-08	Humeral Component, 8", Regular
32-8105-027-04	Humeral Component, 4", Extra-Small
32-8105-027-06	Humeral Component, 6", Extra-Small
32-8105-035-06	Humeral Component, Long Flange, 6", Small
32-8105-035-08	Humeral Component, Long Flange, 8", Small
32-8105-036-06	Humeral Component, Long Flange, 6", Regular
32-8105-036-08	Humeral Component, Long Flange, 8", Regular

Ulnar Implants

Prod. No.	Description
32-8105-053-01	Ulna Component, Plasma, Small, Left
32-8105-053-02	Ulna Component, Plasma, Small, Right
32-8105-063-01	Ulna Component, Plasma, Regular, Left
32-8105-063-02	Ulna Component, Plasma, Regular, Right
32-8105-073-01	Ulna Component, Plasma Long, Small, Left
32-8105-073-02	Ulna Component, Plasma Long, Small, Right
32-8105-083-01	Ulna Component, Plasma Long, Regular, Left
32-8105-083-02	Ulna Component, Plasma Long, Regular, Right
32-8105-093-01	Ulna Component, Plasma Long, Extra-Small, Left
32-8105-093-02	Ulna Component, Plasma Long, Extra-Small, Right
32-8105-043-01	Ulna Component, Plasma, 3-Inch Extra-Small, Left
32-8105-043-02	Ulna Component, Plasma, 3-Inch Extra-Small, Right

Replacement Kits* for Interchangeable System

Prod. No.	Description
32-8105-025-01	Pin/Bushing Replacement Kit Small/Regular
32-8105-027-01	Pin/Bushing Replacement Kit Extra-Small

Revision Old/New Interchange Kits*

Prod. No.	Description
32-8105-029-00	Old Regular Humeral/New Ulna Revision Kit
32-8105-030-00	Old Regular Ulna/New Humeral Revision Kit

Current Replacement/Revision Kits*

Prod. No.	Description
32-8106-000-12	Old Small Bushing Replacement/Revision Kit
32-8106-000-13	Old Regular Bushing Replacement/Revision Kit

* Consult Revision/Replacement Matrix 97-8106-004-00 for proper kit.

Provisionals

Humeral Provisionals

Prod. No.	Description
31-8105-025-04	Humeral, 4", Small
31-8105-025-06	Humeral, 6", Small
31-8105-025-08	Humeral, 8", Small
31-8105-026-04	Humeral, 4", Regular
31-8105-026-06	Humeral, 6", Regular
31-8105-026-08	Humeral, 8", Regular
31-8105-027-04	Humeral, 4", Extra-Small
31-8105-027-06	Humeral, 6", Extra-Small
31-8105-035-06	Humeral, 6", Long Flange, Small
31-8105-035-08	Humeral, 8", Long Flange, Small
31-8105-036-06	Humeral, 6", Long Flange, Regular
31-8105-036-08	Humeral, 8", Long Flange, Regular

Ulnar Provisionals

Prod. No.	Description
31-8105-053-01	Ulna, Small, Left
31-8105-053-02	Ulna, Small, Right
31-8105-063-01	Ulna, Regular, Left
31-8105-063-02	Ulna, Regular, Right
31-8105-073-01	Ulna, Long, Small, Left
31-8105-073-02	Ulna, Long, Small, Right
31-8105-083-01	Ulna, Long, Regular, Left
31-8105-083-02	Ulna, Long, Regular, Right
31-8105-093-01	Ulna, Long, Extra-Small, Left
31-8105-093-02	Ulna, Long, Extra-Small, Right
31-8105-043-01	Ulna, 3-Inch Extra-Small, Left
31-8105-043-02	Ulna, 3-Inch Extra-Small, Right

1. Data on file with Data Monitor, New York, NY 10016-5802 U.S.A.
2. Morrey BF, Adams RA, Bryan RS. Total replacement for post-traumatic arthritis of the elbow. *J Bone Joint Surg (Br)*. 1991;73-B(4):607-612.
3. Morrey BF, Adams RA. Semiconstrained arthroplasty for the treatment of rheumatoid arthritis of the elbow. *J Bone Joint Surg*. 1992;74-A(4):479-490.
4. Morrey BF, Adams RA. Semiconstrained elbow replacement for distal humeral nonunion. *J Bone Joint Surg (Br)*. 1995;77-B:67-72.
5. O'Driscoll SW. Prosthetic elbow replacement for distal humeral fractures and nonunions. *Op Tech in Orthop*. 1994;4(1):54-57.
6. Gill DR, Morrey BF. The Coonrad/Morrey total elbow arthroplasty in patients who have rheumatoid arthritis. *J Bone Joint Surg*. 1998;89(9):1327-1335.

Instruments

Prod. No.	Description
31-8106-027-00	Humeral Guide, Extra-Small
31-8106-040-00	Hinge Pin Removal Tool
31-8106-051-00	Humeral Yolk Gauge
31-8106-052-00	Elbow Extractor
31-8106-053-01	Ulnar Rasp, Small, Left
31-8106-053-02	Ulnar Rasp, Small, Right
31-8106-063-01	Ulnar Rasp, Regular, Left
31-8106-063-02	Ulnar Rasp, Regular, Right
31-8106-065-00	Humeral Rasp, Regular
31-8106-066-00	Humeral Rasp, Small
31-8106-067-00	Starter Ulnar Rasp, Small, Left
31-8106-068-00	Starter Ulnar Rasp, Small, Right
31-8106-069-00	Pilot Ulnar Rasp
31-8106-070-00	Elbow Impactor
31-8106-076-00	Humeral Rasp, 4", Small
31-8106-077-00	Humeral Rasp, 4", Regular
31-8106-078-00	Humeral Starter Rasp
31-8106-079-00	Humeral Starter Rasp, 4"
31-8106-080-00	T-Handle
31-8106-081-00	Humeral Alignment Guide
31-8106-095-00	Humeral Cutting Guide Small/Regular
31-8106-167-00	Starter Awl, 2mm diameter
31-8106-168-00	Starter Awl, 3mm diameter
31-8106-169-00	Starter Awl, 5mm diameter
00-6601-036-00	Awl Reamer
00-8106-090-00	Elbow Instrument Case

SETS

Prod. No.	Description
31-8106-000-01	Coonrad/Morrey Instrument Set & Case
31-8106-000-02	Long Ulnar Provisional Set
31-8106-000-03	Long Flange Provisional Set
31-8106-000-04	Extra-Small Humeral Instrument Set
31-8106-000-06	Extra-Small Ulna Provisional Set
32-8105-000-00	Standard Implant Set
32-8105-000-02	Long Ulnar Implant Set
32-8105-000-03	Long Flange Humeral Implant Set
32-8105-000-04	Extra-Small Humeral Implant Set
32-8105-000-05	Bushing Replacement Kit
32-8105-000-06	Extra-Small Ulna Implant Kit Set



Contact your Zimmer representative or visit us at www.zimmer.com

